


Dr Mark's HyGenie®

7 Simple Steps For Peak Oral Hygiene


1 OPEN

Hold the bottom handle (blue) and device body (clear) in one hand. Turn the top handle (white) either direction until it clicks into place. Lift off the top.

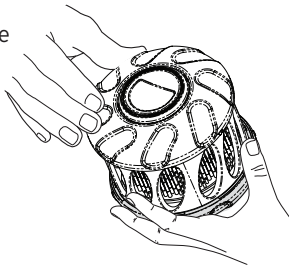


3 CLOSE

Line up the arrows found on the top handle and clear body, as shown. Gently press down until the handle 'sits' into place. Twist it in either direction and it will 'click & lock' closed.

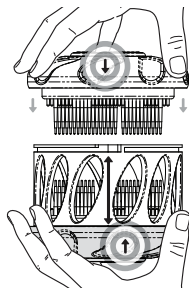

5 WASH


Rotate the top and bottom handles in opposite directions, or, alternate top and bottom washing by holding one handle in place and rotating the opposite handle. Wash for 30-45 seconds or until DentalFresh foam thoroughly covers your appliance.

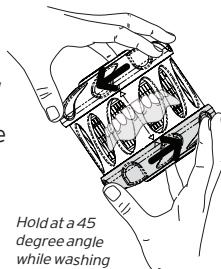

Hold at a 45 degree angle while washing as illustrated


7 WEAR IT, STORE IT, PROTECT IT

Now that your dental appliance is clean, fresh and ready, you can wear it or store it safely in Dr Mark's HyGenie® until you need it next.

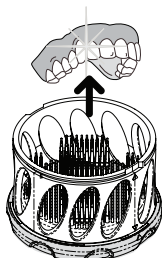


2 PLACEMENT

Nestle your dental appliance into the bottom brush. Paired orthodontic aligners can be placed together, with one facing upwards and one facing downwards, as shown.

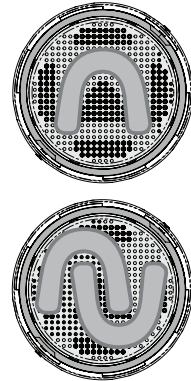


4 DENTAL FRESH

Holding your HyGenie over a basin, squeeze 1 to 2 generous drops of Dr Mark's DentalFresh® into the device through the perforations and apply a splash of water.

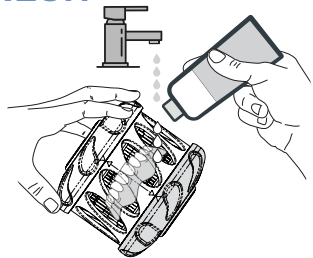

6 RINSE

Rinse the HyGenie under gently running water, rotating the handles in opposite directions, until all DentalFresh foam is rinsed away and water runs clear through the device.

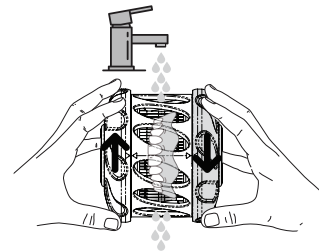

For more Information:



Whip Mix
800-626-5651
www.whipmix.com